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www.respirox.com.tr

CE0123

The model of the Heated Humidifier is BMC-50. The Heated Humidifier is designed only for use with specific Ece Respirox® CPAP / Auto CPAP or BPAP devices. Do not use with any other devices.

The humidifier moistens the air delivered by the Ece Respirox® CPAP / BPAP devices.

The Heated Humidifier is only used for single patient and must not be re-used on another person. This is to avoid the risk of cross-infection.

The Heated Humidifier is not intended for use with a patient whose upper airways have been bypassed.

Table of Contents

Intended Use	1
Warnings & Cautions.....	2
Symbols.....	3
Features	4
Set Up	4
Daily Use	5
Cleaning	5
Disinfection	6
Service	6
Specifications	7
Disposal	7
Traveling with the System.....	7
EMC Requirements.....	8
Warranty.....	13

Intended Use

Heated Humidifier is designed for specific ECE CPAP / Auto CPAP or BPAP devices and is used to increase the outlet air humidity. These devices are intended to deliver positive pressure for the treatment of adult Obstructive Sleep Apnea (OSA) only, either in the hospital or at home.

Contraindications

1. When Heated Humidifier is used along with the Respirox® CPAP / Auto CPAP System, the contraindications are as follows:

Studies have shown that the following pre-existing conditions may contraindicate the use of positive airway pressure therapy for some patients:

- Bullous Lung Disease
- Bypassed Upper Airway
- Pneumothorax
- Pathologically Low Blood Pressure
- Pneumocephalus has been reported in a patient using nasal Continuous Positive Airway Pressure. Caution should be used when prescribing CPAP for susceptible patients such as those with: cerebral spinal fluid leaks, abnormalities of the cribriform plate, prior history of head trauma, and / or pneumocephalus.

(Chest 1989; 96:1425-1426)

- The use of positive airway pressure therapy may be temporarily contraindicated if you exhibit signs of a sinus or middle ear infection.

Contact your health care professional if you have any questions concerning your therapy.

2. When Heated Humidifier is used along with the Respirox® BPAP System, the contraindications are as follows:

The BPAP should not be used if you have an insufficient respiratory drive to endure brief interruptions in non-invasive ventilation therapy. The BPAP is not a life support ventilator and may stop operating with power failure or in the unlikely event of certain fault conditions.

If you have any of the following conditions, tell your doctor before using the BPAP:

- acute sinusitis or otitis media.
- conditions predisposing to a risk of aspiration of gastric contents.
- epistaxis causing a risk of pulmonary aspiration.
- hypotension or significant intravascular volume depletion.
- Inability to maintain a patent airway or adequately clear secretions.
- pneumothorax or pneumomediastinum.
- recent cranial trauma or surgery.

Warnings & Cautions

IMPORTANT!

- Read all instructions before using the humidifier.
- Use only with Ece Respirox® devices whose instructions specify the use of this humidifier.
- Please use the mask which meets ISO 17510:2015.

CAUTIONS!

- Indicates the possibility of damage to the device.
- US federal law restricts this device to sale by or on the order of a physician.
- If fluids are spilled onto the humidifier platform, unplug the power cord from the AC wall outlet and allow the humidifier platform to drain and dry before using.
- Take precautions to protect furniture from water damage.

WARNINGS!

- Indicates the possibility for injury to the user or the operator.
- Use the humidifier only for its intended use as described in this manual.
- Use only accessories recommended by Ece.
- Never operate the humidifier if any of the parts are damaged, if it is not working properly, or if the humidifier has been dropped or mishandled. Do not use the humidifier if the water chamber is leaking or damaged in any way. Have any damaged parts replaced before continuing use.
- Never touch the heater plate unless the humidifier is unplugged and the plate has cooled down.
- This equipment is not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide.
- Periodically inspect the power cord for signs of wear or damage. Replace if necessary.
- Do not operate the device in direct sunlight or near a heating appliance because these conditions can increase the temperature of the air coming out of the device.
- When humidifier is used outside the specified ambient temperature range or humidity range, the performance of humidifier will be compromised.
- Do not bring the device or accessories into a Magnetic Resonance (MR) environment as it may cause unacceptable risk to the patient or damage to the device or MR medical devices. The device and accessories have not been evaluated for safety in an MR environment.
- Do not use the device or accessories in an environment with electromagnetic equipment such as CT scanners, Diathermy, RFID and electromagnetic security systems (metal detectors) as it may cause unacceptable risk to the patient or damage to the device. Some electromagnetic sources may not be apparent, if you notice any unexplained changes in the performance of this device, if it is making unusual or harsh sounds, disconnect the power cord and discontinue use. Contact your home care provider.

Symbols

Follow Instructions for Use

Type BF Applied Part (mask)

Class II (Double Insulated)

DC Power

IP21

≥ 12.5 mm Diameter, Drip-Proof, Vertical

Serial Number of the Product

Manufacturer


Authorized Representative in the European Community

European CE Declaration of Conformity

Features

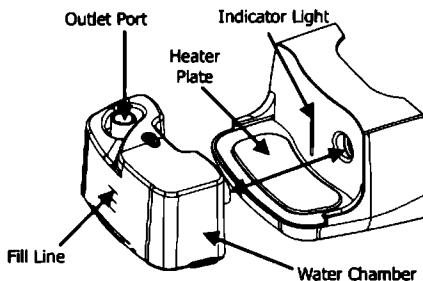
Fill Line: This indicates the maximum water level for safe operation.

Outlet Port: Connect the flexible tubing (coaxial 22 mm) here.

Water Chamber: The removable water chamber holds the water for humidification. It has a silica gel cover to outflow the water in cleaning.

Heater Plate: Warms the water in the water chamber.

Indicator Light: When lit, this indicates that the heater plate has been turned on.

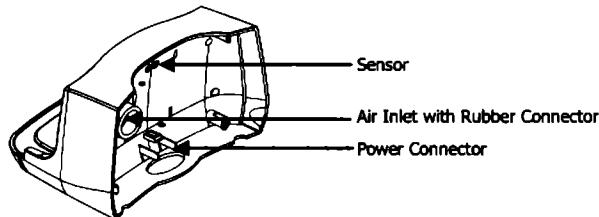
Sensor: Do not block this sensor. The humidifier will not operate unless it can detect that the Respirox® device is connected correctly.

Air Inlet with Rubber Connector: Connect to the outlet port on the Respirox® device.

Power Connector: Plug this connector into the power outlet on the Ece Respirox® device.

Set Up

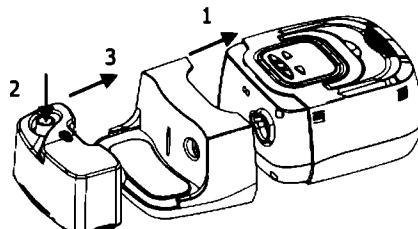1. Connect the humidifier to the Respirox® device. Make sure they are connected completely.
2. Place the Respirox® device and humidifier on a firm, flat surface at a level lower than your sleeping position.

CAUTIONS!

- Do not turn the humidifier on without the water chamber installed.
- Do not touch the metal parts of the water chamber during device operation, otherwise it may cause burns.
- Take precautions to protect furniture from water damage.

Daily Use

1. Fill the chamber to the fill line with water (approx. 350 ml) from the outlet port. Distilled water is recommended. Do not overfill the water chamber.

CAUTION!

- Always remove the chamber from the humidifier before filling with water.
2. Press down the spring loaded heater plate with the water chamber and slide the chamber into place. Make sure the rubber connector on the inlet port fits securely over the Respirox® device's air outlet.
3. Connect the flexible tubing to the outlet port on the water chamber.

CAUTION!

- Avoid moving or tilting the humidifier when the water chamber has water in it.
4. When the Respirox® device begins blowing, the humidifier will work automatically. The yellow indicator light on the humidifier will turn on. Press the humidifier button  can turn off or restart the humidifier.
5. The ideal humidity setting depends on room temperature and humidity. Initially, a setting of 3 is recommended. You can adjust this setting at any time. Please adjust the humidifier setting according to the CPAP User Manual.

IMPORTANT!

- When the air flow is turned off, the humidifier will automatically shut off.

Cleaning

The water chamber should be cleaned daily or after each use.

WARNING!

- Emptying and cleaning the water chamber daily will help to prevent mold and bacteria growth.

WARNING!

- Allow the water in the chamber to cool down to room temperature before removing it from the humidifier.

WARNING!

- To avoid electrical shock, disconnect the power cord of Respirox® device before cleaning the humidifier. DO NOT immerse the humidifier into any fluids.

Water Chamber

1. Turn the Respirox® device off and allow approximately 15 minutes for the heater plate and water to cool.
2. Disconnect the tubing from the water chamber. Press down on the water chamber and slide it out of the humidifier platform.
3. Open the silica gel cover of water chamber and discard any remaining water. Fill a solution of

warm water and a mild dishwashing detergent into the chamber, cover the silica gel cover, rock the chamber a few minutes, and then outflow the solution. Rinse the chamber several times with clean water and allow to air dry.

4. Fill the water chamber and close the silica gel cover. Inspect the water chamber for any leak or damage. Replace the water chamber if any damage is present.

Humidifier Platform

1. Clean the humidifier platform by wiping with a damp cloth. Allow to air dry.
2. Inspect the humidifier platform for any damage and replace if necessary.

Silicone Tubing of Humidifier

Clean the Silicone Tubing by wiping with an alcohol cotton stick. Allow to air dry.

Disinfection

Generally speaking, if you have strictly followed the above cleaning instructions, you do not have to disinfect the humidifier. If the device is contaminated or used in clinical trials, you may purchase disinfectants from a pharmacist to disinfect the water chamber.

Disinfection of Humidifier Water Chamber

Prior to disinfection, clean the water chamber according to Section "Cleaning". The disinfection methods are as follows:

- (1) Heat disinfection: Disinfect the water chamber by immersing it in tap water at $75^{\circ}\text{C} \pm 2^{\circ}\text{C}$ for 30 minutes.
- (2) Use mild disinfectants.

CAUTIONS!

- Disinfectants tend to damage materials and reduce the life of components. Try to select the appropriate disinfectant, and follow the disinfectant manufacturer's instructions and recommendations.
- After disinfection, check the disinfected component for any signs of damage. Replace any damaged component immediately.

WARNINGS!

- After disinfection, rinse any disinfected component in clean water thoroughly, especially components in close contact with the patient such as the mask, headgear, and tube, so as to prevent disinfectant residuals from damaging the skin or respiratory tract or causing allergies.
- The device shall not be serviced or maintained while in use with a patient.
- Sterilization of this device and its components other than recommended is not permitted.

Service

The humidifier does not require routine servicing.

If the humidifier malfunctions, contact your home care provider immediately. Never attempt to open the humidifier's enclosure. If necessary, contact your local authorized dealer or ECE TİBBİ CİHAZLAR VE MEDİKAL SAN.TİC. A.Ş, for technical support and documents.

Specifications

Size

Dimensions: 120 mm x 194 mm x 112 mm

Weight: 0.8 kg

Water Capacity: 350 ml at recommended water level

Product Use, Transport and Storage

	Operation	Transport and Storage
Temperature:	5°C to 35°C	-20°C to 55°C
Humidity:	≤ 80% Non-condensing	≤ 93% Non-condensing
Atmospheric Pressure:	700 to 1060 hPa	500 to 1060 hPa

Power Requirements

24 V DC 1.0 A max

Type of Protection Against Electric Shock

Class II Equipment

Degree of Protection Against Electric Shock

Type BF Applied Part

Degree of Protection Against Ingress of Water

IP21 – ≥ 12.5 mm Diameter, Drip-Proof, Vertical

Heater Settings

1 to 5 (104°F to 149°F / 40°C to 65°C)

Maximum Operating Pressure

40 hPa

Pressure Drop with Humidifier

< 0.5 hPa at 60 LPM flow

Maximum Delivered Gas Temperature

< 40°C

Humidity Range

10 to 45 mg H₂O/L

Disposal

When necessary, dispose of the device and accessories in accordance with local regulations.

Traveling with the System

Packing the System

When traveling, the optional CPAP carrying case is for carry-on luggage only. The carrying case will not protect the humidifier if it is put through checked baggage.

Security Stations

For ease at security stations, there is a note on the bottom of the humidifier stating that it is medical equipment. It may be helpful to bring this manual along with you for security personnel.

EMC Requirements

Guidance and manufacturer's declaration - electromagnetic emissions		
The device is intended for use in the electromagnetic environment specified below. The user of the device should ensure that it is used in such an environment.		
Emissions Test	Compliance	Electromagnetic Environment - Guidance
RF emissions CISPR 11	Group 1	The device uses RF energy only for its internal function. Therefore its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Class A	The device is suitable for use in all establishments including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes
Voltage fluctuations / flicker emissions IEC 61000-3-3	Complies	

Guidance and manufacturer's declaration - electromagnetic immunity			
The device is intended for use in the electromagnetic environment specified below. The user of the device should make sure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 8 kV contact ± 15 kV air	± 8 kV contact ± 15 kV air	Floor should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%
Electrical fast transient / burst IEC 61000-4-4	± 2 kV for power supply lines	± 2 kV for power supply lines	Mains power quality should be that of a typical commercial or hospital environment
Surge IEC 61000-4-5	± 1 kV line (s) to line (s)	± 1 kV line (s) to line (s)	Mains power quality should be that of a typical commercial or hospital environment
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% U_T , 0.5 cycle At 0° , 45° , 90° , 135° , 180° , 225° , 270° and 315° 0% U_T , 1 cycle 70% U_T , 25 / 30 cycle At 0° 0% U_T , 250 / 300 cycle	0% U_T , 0.5 cycle At 0° , 45° , 90° , 135° , 180° , 225° , 270° and 315° 0% U_T , 1 cycle 70% U_T , 25 / 30 cycle At 0° 0% U_T , 250 / 300 cycle	Mains power quality should be that of a typical commercial or hospital environment. If the user of the device requires continued operation during power mains interruptions, it is recommended that the device be powered from an uninterruptible power supply or a battery
Power frequency (50 / 60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment

Note: U_T is the a.c. mains voltage prior to application of the test level.

Guidance and manufacturer's declaration - electromagnetic immunity			
The device is intended for use in the electromagnetic environment specified below. The user of the device should make sure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Conducted RF IEC 61000-4-6	3 V 0.15 MHz ~ 80 MHz 6 V in ISM and amateur radio bands between 0.15 MHz and 80 MHz	3 V 0.15 MHz ~ 80 MHz 6 V in ISM and amateur radio bands between 0.15 MHz and 80 MHz	Portable and mobile RF communications equipment should be used no closer to any part of the device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.7 GHz	10 V/m 80 MHz to 2.7 GHz	<p>Recommended separation distance</p> $d = 1.17\sqrt{P}$ $d = 0.35\sqrt{P} \quad 80 \text{ MHz to } 800 \text{ MHz}$ $d = 0.70\sqrt{P} \quad 800 \text{ MHz to } 2.5 \text{ GHz}$ <p>Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,^a should be less than the compliance level in each frequency range.^b</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
<p>Note 1: At 80 MHz and 800 MHz, the higher frequency range applied.</p> <p>Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.</p>			
<p>^a Field strengths from fixed transmitters, such as base stations for radio (cellular / cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the device is used exceeds the applicable RF compliance level above, the device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the device.</p> <p>^b Over the frequency range 150 kHz to 80 MHz, the field strengths should be less than 10 V/m.</p>			

Recommended separation distances between portable and mobile RF communications equipment and the device

The device is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the device can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the device as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output of transmitter W	150 kHz ~ 80 MHz $d = 1.17\sqrt{P}$	80 MHz ~ 800 MHz $d = 0.35\sqrt{P}$	800 MHz ~ 2.5 GHz $d = 0.70\sqrt{P}$
0.01	0.12	0.04	0.07
0.1	0.37	0.12	0.23
1	1.17	0.35	0.70
10	3.70	1.11	2.22
100	11.7	3.50	7.00

Note 1: At 80 MHz and 800 MHz, the higher frequency range applied.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Recommended separation distances between RF wireless communications equipment						
Frequency MHz	Maximum Power W	Distance	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance	
385	1.8	0.3	27	27	RF wireless communications equipment should be used no closer to any part of the device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $E = \frac{6}{d} \sqrt{P}$ Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitter, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range. Interference may occur in the vicinity of equipment marked with the following symbol: 	
450	2	0.3	28	28		
710	0.2	0.3	9	9		
745						
780	2	0.3	28	28		
810						
870	2	0.3	28	28		
930						
1720	2	0.3	28	28		
1845						
1970	2	0.3	28	28		
2450						
5240	0.2	0.3	9	9		
5500						
5785						

Note: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

WARNINGS!

- This device should not be used in the vicinity or on the top of other electronic equipment such as cell phone, transceiver or radio control products. If you have to do so, the device should be observed to verify normal operation.
- The use of accessories and power cord other than those specified, with the exception of cables sold by the manufacturer of the equipment or system as replacement parts for internal components, may result in increased emissions or decreased immunity of the equipment or system.

Warranty

ECE TİBBİ CİHAZLAR VE MEDİKAL SAN.TİC. A.Ş warrants that this humidifier shall be free from defects of workmanship and materials and will perform in accordance with the product specifications for a period of one (1) year from the date of sale by ECE TİBBİ CİHAZLAR VE MEDİKAL SAN.TİC. A.Ş to the dealer. If the product fails to perform in accordance with the product specifications, ECE TİBBİ CİHAZLAR VE MEDİKAL SAN.TİC. A.Ş will repair or replace, at its option, the defective material or part. ECE TİBBİ CİHAZLAR VE MEDİKAL SAN.TİC. A.Ş will pay customary freight charges from ECE TİBBİ CİHAZLAR VE MEDİKAL SAN.TİC. A.Ş to the dealer location only. This warranty does not cover damage caused by accident, misuse, abuse, alteration and other defects not related to material or workmanship.

To exercise your rights under this warranty, contact your local, authorized dealers or:

MANUFACTURER:

BMC Medical Co., Ltd.

Room 110 Tower A Fengyu Building, No. 115 Fucheng Road, Haidian, 100036
Beijing, PEOPLE'S REPUBLIC OF CHINA

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Tel: +90 530 306 44 39

Faks: +90 212 528 07 17

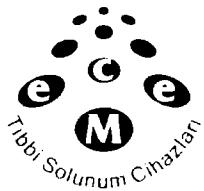
eposta volkancalman@respirox.com.tr

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